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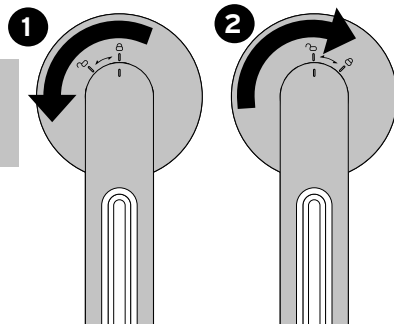
en

Safety warnings

Read the safety warnings carefully and only use the product as described in these instructions to avoid accidental injury or damage.

Intended use

- The product is designed for removing hard skin that forms on your feet. It can also be used for removing hard skin that forms on hands and elbows.
- The product is designed for domestic use and is not suitable for commercial purposes, in particular not for medical foot care.
- Do not use the device...

... if you suffer from diabetes or other wound-healing disorders. The possibility of injuring yourself with this device is very slight. Nevertheless, we are obliged to draw your attention to the fact that you should not use it if you suffer from the illnesses mentioned above.

... on reddened or damaged skin or on burns. Moreover, do not use the device on areas of the body without hard skin, on birthmarks, moles or warts, or on areas affected by skin disorders.

... if another person has already used it. The device and the accompanying attachments should be used by one person only. This will prevent the spread of infections such as athlete's foot.

DANGER to children

- Keep children away from the packaging material. It poses several risks, including the risk of suffocation!
- This device must not be used on infants due to the sensitivity of their skin.

WARNING - risk of injury

- Stop removing hard skin immediately if you start to feel pain!
- Do not rub the device on one area for a prolonged period; the resulting friction will otherwise generate intense heat in that area.

CAUTION - material damage

- Stop using the attachment if it is very worn down.
- Clean the attachments thoroughly after every use as outlined in the chapter “Cleaning and disinfecting the device”. This is to remove grease and other residues.
- Do not use any caustic, aggressive or abrasive products for cleaning.

Please note

Hard skin is normal, especially on feet. However, it should not be allowed to become too thick or hard as this could cause pain. Excessive formation of hard skin can be caused by shoes that are too small or do not fit properly, by foot deformities, or by extremely dry skin. If you suffer from excessive hard skin, try to treat the causes and consult a doctor for advice.

To prevent the excessive formation of hard skin, we recommend that you take regular foot baths and then apply lotion or foot creams to your feet. Make sure that your shoes fit well and try to wear different shoes every day if possible.

Use

- We recommend beginning with the pumice attachment to remove dead skin, then using the callus remover attachment, and finally using the brush attachment to clean your feet. You can of course also use the attachments individually.
- Remove too little skin rather than too much. You should always leave a small amount of hard skin on your feet to ensure that your feet remain protected. We recommend that you remove hard skin regularly but gently.

Changing attachments

1. Turn the handle towards then detach it.
2. Attach the handle to the required attachment, ensuring that the line marking on the handle is aligned with the line marking of the symbol on the attachment, then turn the handle towards until you feel and hear it click into place.

Removing hard skin

1. To achieve better results, take a 5- to 10-minute foot bath beforehand. Always dry your feet thoroughly after bathing.
2. Use your hand to tauten the area of the skin you wish to treat.
3. Applying gentle pressure, carefully move the attachment back and forth over the hard skin you wish to remove.
4. Brush the removed skin off using the brush attachment.
5. Rub a suitable body cream into the area of skin treated. Do not use care products that might irritate your skin.
6. After use, clean the attachments as described in the chapter “Cleaning and disinfecting the device”.

Cleaning and disinfecting the device

- ▷ Brush the rough surfaces of the pumice and callus remover attachments with a small soft brush.
- ▷ If necessary, wipe the other parts with a slightly dampened cloth.
- ▷ If necessary, disinfect the attachments using ethyl alcohol, which is available from chemists.
- ▷ Allow the product to dry thoroughly before putting it away.

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